



How to **transfer** experiments

from

Signals Notebook **Individual Edition**

to

Signals Notebook **Standard**




As an Administrator


If you are not a Signals Notebook Admin, go directly to page 3.

Note that you **CANNOT** import your Individual Edition Notebooks into Signals Notebook Standard **UNLESS** your Signals Notebook Standard Admin has gone through these steps.

Step 1:
In your user profile, go to **System Configuration**

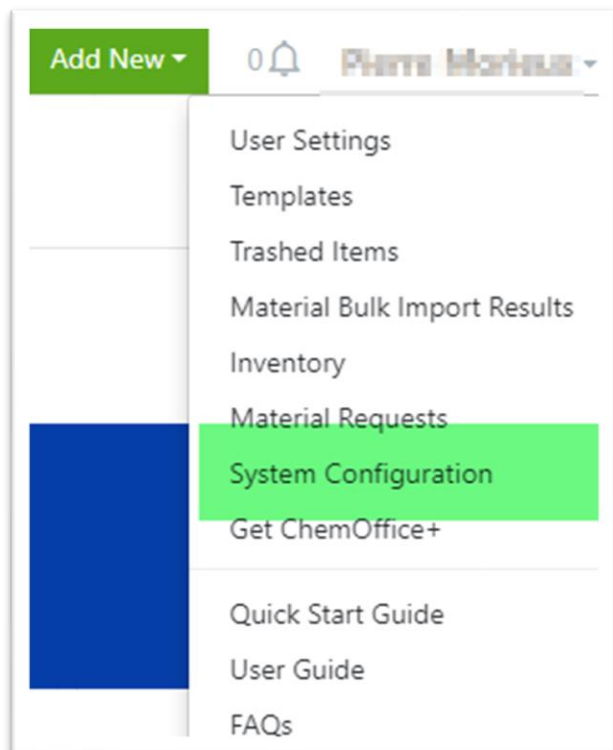


Step 2:
In the System Configuration go to **System Settings**

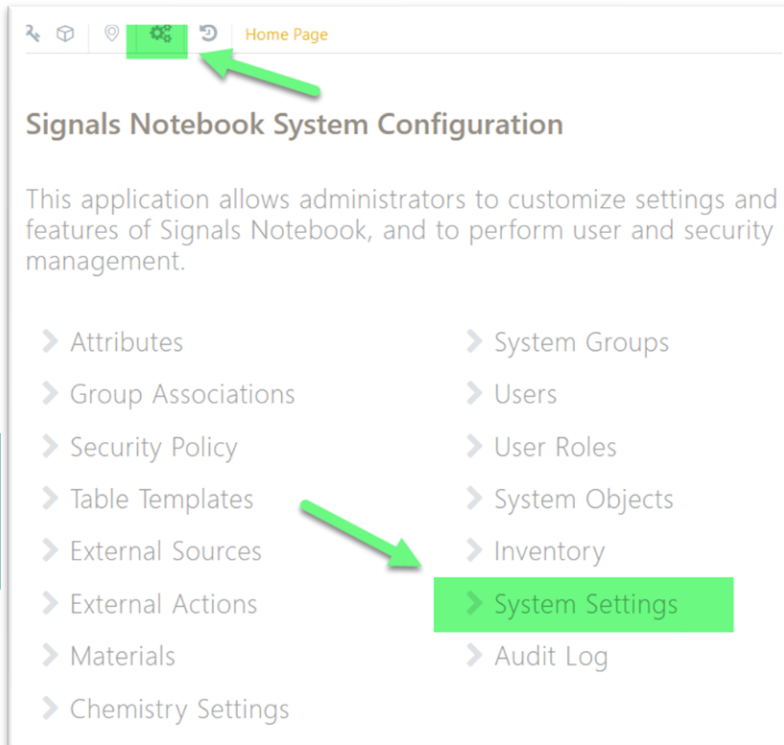

Step 3:
In System Settings > Other > **Workflow Settings**

Toggle the « **Enable Import of Notebooks from Signals Notebook Individual** » Option

Save the Settings

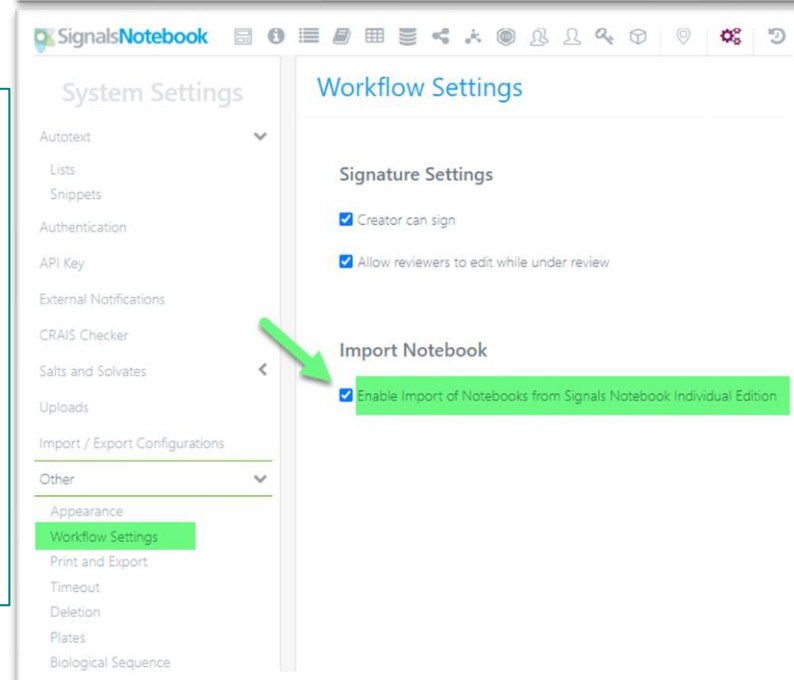



As a Signals Notebook **Individual Edition** User


- You can only export **NOTEBOOKS** and their content from Signals Notebook **Individual Edition**.
- Any **EXPERIMENT** outside of a **NOTEBOOK** needs to be moved into a **NOTEBOOK** for export.
- Exported **NOTEBOOKS** can only contain **up to 100 EXPERIMENTS**. If your **NOTEBOOKS** have **more than 100 EXPERIMENTS**, they will need to be broken down during the export process (Notebook-0-to-100, Notebook-101-to-200, Notebook-201-to-300 etc.)
- **EXPERIMENT** hyperlinks, user comments, and audit trail information from Signals Notebook **Individual Edition** will not be exported.
- You can only move **NOTEBOOKS** from <https://chemdraweln.signalsresearch.revivitycloud.com> towards a **SNB Standard** instance or from <https://signalsresearch.revivitycloud.com/> towards a **SNB Standard** instance
- You **CANNOT move data** between <https://chemdraweln.signalsresearch.revivitycloud.com> and <https://signalsresearch.revivitycloud.com/>



As a Signals Notebook User


Log on to Signals Notebook **Individual Edition**

Go into a **NOTEBOOK** you wish to export.
In the 3 dots menu, select « Export to File »

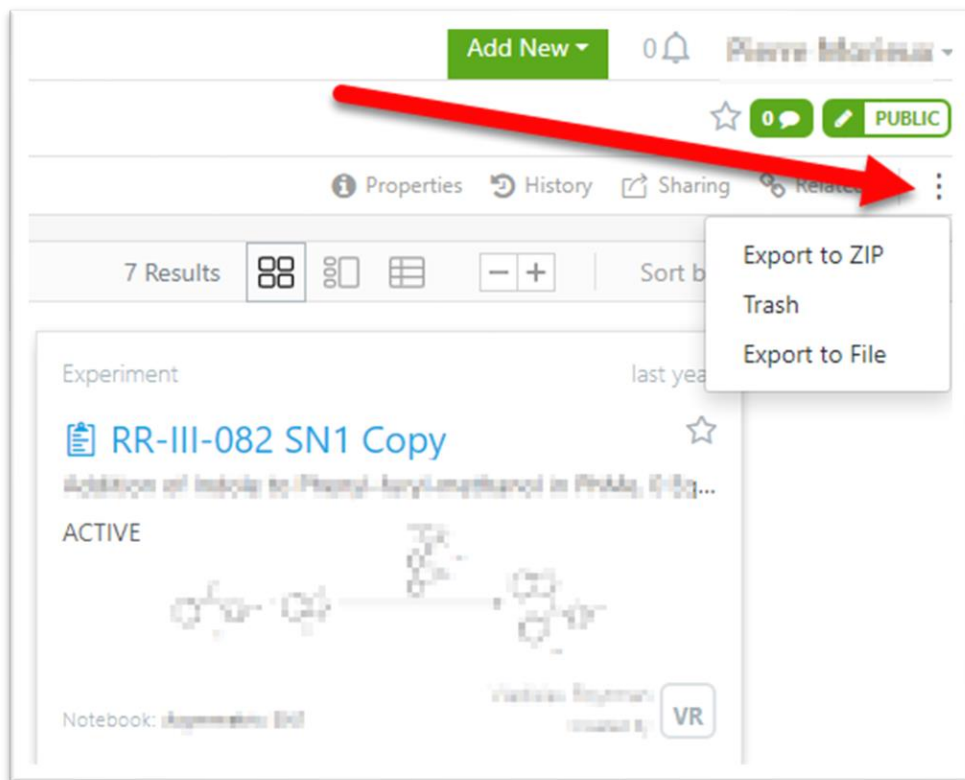

Your **NOTEBOOK** will get exported as a **.snb.notebook** file

Repeat for every **NOTEBOOK**

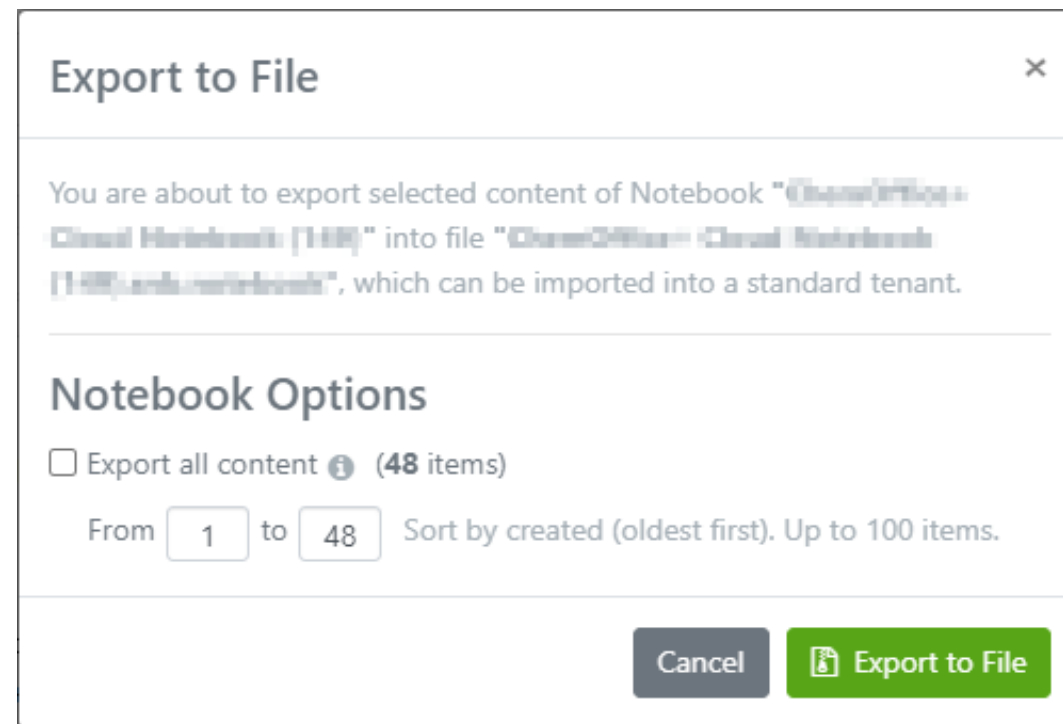


Note: if your Notebook contains more than 100 experiments, you will need to export it into multiple files.



As a Signals Notebook User

Log on to Signals Notebook **Standard**

Under your User Profile, click on
« Import Notebook from File »

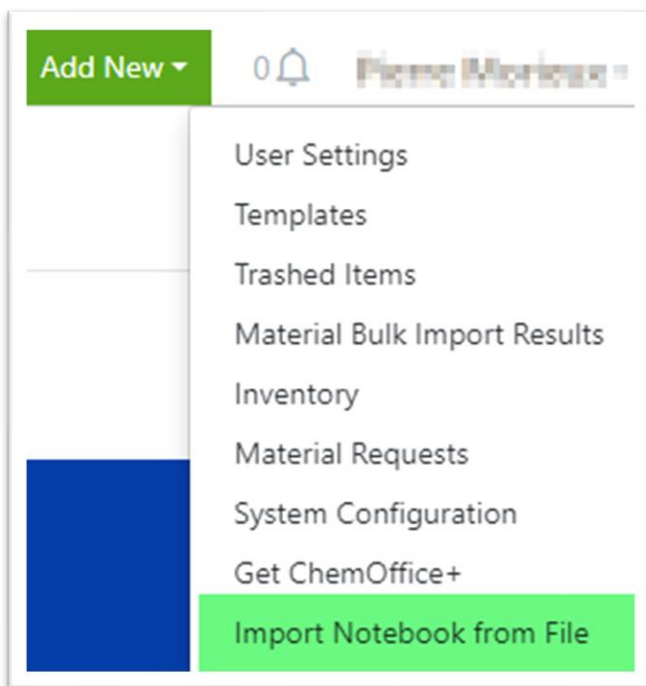

If you cannot see this option, reach out to your Admin

Select the *.snb.notebook* file that you wish to upload
to Signals Notebook **Standard**

Click on **Import**

